

Evaluasi biokompatibilitas saluran gas pernapasan pada penerapan pelayanan kesehatan - Bagian 2: Uji emisi partikulat

(ISO 18562-2:2017, IDT, Eng)

© ISO 2017 – All rights reserved

© BSN 2020 untuk kepentingan adopsi standar © ISO menjadi SNI - Semua hak dilindungi

Hak cipta dilindungi undang-undang. Dilarang mengumumkan dan memperbanyak sebagian atau seluruh isi dokumen ini dengan cara dan dalam bentuk apapun serta dilarang mendistribusikan dokumen ini baik secara elektronik maupun tercetak tanpa izin tertulis dari BSN

BSN

Email: dokinfo@bsn.go.id

www.bsn.go.id

Diterbitkan di Jakarta

Daftar isi

Daftar isi	i
Prakata	ii
Introduction	iii
1 Scope	1
2 Normative reference	2
3 Term and definition	2
4 General principles	2
5 * PARTICULATE MATTER emissions	3
Annex A (informative) Rationale and guidance	8
Annex B (informative) Reference to the essential principles	13
Annex C (informative) Terminology — Alphabetized index of defined terms	14
Bibliography	15

Prakata

Standar Nasional Indonesia (SNI) ISO 18562-2:2017, dengan judul *Evaluasi biokompatibilitas saluran gas pernapasan pada penerapan pelayanan kesehatan - Bagian 2: Uji emisi partikulat (ISO 18562-2:2017, IDT, Eng)*, merupakan hasil adopsi identik dari standar ISO 18562-2:2017 *Biocompatibility evaluation of breathing gas pathways in healthcare applications — Part 2: Tests for emissions of particulate matter*, dengan metode republikasi *reprint*, yang ditetapkan oleh BSN pada tahun 2020.

Standar ini disusun oleh Komite Teknis 11-03 Alat Kesehatan Elektromedik dengan Badan Standardisasi Nasional (BSN) sebagai sekretariat Komite Teknis. Standar ini telah dibahas dalam rapat teknis, dan terakhir disepakati dalam rapat konsensus di Jakarta pada tanggal 20 April 2020 yang dihadiri oleh para pemangku kepentingan (*stakeholder*) terkait, yaitu perwakilan dari produsen, konsumen, pakar dan pemerintah, serta perwakilan dari lembaga penguji, asosiasi, perguruan tinggi, pakar serta instansi terkait.

Standar ini telah melalui tahap jajak pendapat pada tanggal 18 Mei 2020 sampai dengan 6 Juni 2020 dengan hasil akhir disetujui menjadi SNI.

Apabila di kemudian hari pengguna menemukan kesulitan dalam penggunaan standar ini, maka dianjurkan untuk merujuk pada standar aslinya yaitu ISO 18562-2:2017 dan/atau dokumen terkait lain yang menyertainya.

Perlu diperhatikan bahwa kemungkinan beberapa unsur dari dokumen standar ini dapat berupa hak paten. Badan Standardisasi Nasional tidak bertanggungjawab untuk pengidentifikasian salah satu atau seluruh hak paten yang ada.

Introduction

This document is intended to protect PATIENTS connected to MEDICAL DEVICES from excessive amounts of PARTICULATE MATTER that arises from within GAS PATHWAYS of MEDICAL DEVICES.

This document is intended to cover the biological evaluation of GAS PATHWAYS of MEDICAL DEVICES within a RISK MANAGEMENT PROCESS, as part of the overall MEDICAL DEVICE evaluation and development. This approach combines the review and evaluation of existing data from all sources with, where necessary, the selection and application of additional tests.

In general, the ISO 10993 series^[2] is intended to cover the biological evaluation of MEDICAL DEVICES. However, the ISO 10993 series does not appropriately address the biological evaluation of the GAS PATHWAYS of MEDICAL DEVICES. For example, the ISO 10993 tests do not evaluate inspired PARTICULATE MATTER.

It is not within the scope of this document to address contamination arising from the source of the breathing gases entering such MEDICAL DEVICES, but rather address only the potential contamination generated from within the MEDICAL DEVICE itself. This contamination might be from the original manufacturing PROCESS or be generated by the MEDICAL DEVICE itself during use.

This document is concerned with PARTICULATE MATTER that could be conveyed to the PATIENT by the breathing gases. The smaller the particle, the deeper into the lungs it can penetrate and the longer it takes the body to eliminate it. Originally, the main health concerns with regard to PARTICULATE MATTER were focused on respiratory health, but now there is emerging evidence of effects on the cardiovascular system as well.

The tests for the presence of PARTICULATE MATTER generated by respiratory MEDICAL DEVICES are based on standard laboratory practice and require no advanced techniques or equipment.

The acceptable levels of contamination are based on worldwide published health data for particulates. It is accepted that there is no point in setting a level that is lower than that found in air that people might breathe every day of their lives.

In this document, the following print types are used:

- requirements and definitions: roman type;
- informative material appearing outside of tables, such as notes, examples and references: in smaller type. Normative text of tables is also in a smaller type;
- *test specifications: italic type;*
- TERMS DEFINED IN CLAUSE 3 OF THIS DOCUMENT OR AS NOTED: SMALL CAPITALS TYPE.

In this document, the conjunctive “or” is used as an “inclusive or” so a statement is true if any combination of the conditions is true.

The verbal forms used in this document conform to usage described in Annex H of the ISO/IEC Directives, Part 2. For the purposes of this document, the auxiliary verb:

- “shall” means that compliance with a requirement or a test is mandatory for compliance with this document;
- “should” means that compliance with a requirement or a test is recommended but is not mandatory for compliance with this document;

- “may” is used to describe a permissible way to achieve compliance with a requirement or test.

An asterisk (*) as the first character of a title or at the beginning of a paragraph or table title indicates that there is guidance or rationale related to that item in Annex A.

The attention of Member Bodies is drawn to the fact that equipment manufacturers and testing organizations may need a transitional period following publication of a new, amended or revised ISO publication in which to make products in accordance with the new requirements and to equip themselves for conducting new or revised tests. It is the recommendation of the committee that the content of this publication be adopted for implementation nationally not earlier than 3 years from the date of publication for equipment newly designed and not earlier than 5 years from the date of publication for equipment already in production.

Evaluasi biokompatibilitas saluran gas pernapasan pada penerapan pelayanan kesehatan - Bagian 2: Uji emisi partikulat

1 Scope

This document specifies tests for the emissions of PARTICULATE MATTER from the GAS PATHWAYS of a MEDICAL DEVICE, its parts or ACCESSORIES, which are intended to provide respiratory care or supply substances via the respiratory tract to a PATIENT in all environments. The tests of this document are intended to quantify particles from 0,2 μm DIAMETER to 10 μm DIAMETER that are emitted by the MEDICAL DEVICE, its parts or ACCESSORIES into the respirable gas stream. This document establishes acceptance criteria for these tests. This document does not address nanoparticles. Insufficient data exist to establish exposure limits for particles less than 0,2 μm in DIAMETER.

NOTE 1 Smaller and larger particles could also present biological HAZARDS, and additional information outside the scope of this document can be needed to meet requirements of some AUTHORITIES HAVING JURISDICTION.

This document therefore adopts the same approach as the US Environmental Protection Agency (EPA) in setting limits based solely on particle size and not their chemistry.

This document addresses potential contamination of the gas stream arising from the GAS PATHWAYS, which is then conducted to the PATIENT.

This document applies over the EXPECTED SERVICE LIFE of the MEDICAL DEVICE in NORMAL USE and takes into account the effects of any intended processing or reprocessing.

This document does not address biological evaluation of the surfaces of GAS PATHWAYS that are in direct contact with the PATIENT. The requirements for direct contact surfaces are found in the ISO 10993 series.

MEDICAL DEVICES, parts or ACCESSORIES, containing GAS PATHWAYS that are addressed by this document, include, but are not limited to, ventilators, anaesthesia workstations (including gas mixers), breathing systems, oxygen conserving devices, oxygen concentrators, nebulizers, low-pressure hose assemblies, humidifiers, heat and moisture exchangers, respiratory gas monitors, respiration monitors, masks, mouth pieces, resuscitators, breathing tubes, breathing systems filters, Y-pieces, and any breathing ACCESSORIES intended to be used with such devices. The enclosed chamber of an incubator, including the mattress, and the inner surface of an oxygen hood are considered to be GAS PATHWAYS and are also addressed by this document.

This document does not address contamination already present in the gas supplied from the gas sources while MEDICAL DEVICES are in NORMAL USE.

EXAMPLE Contamination arriving at the MEDICAL DEVICE from gas sources such as MEDICAL GAS PIPELINE SYSTEMS (including the non-return valves in the pipeline outlets), outlets of pressure regulators connected or integral to a medical gas cylinder, or room air taken into the MEDICAL DEVICE is not addressed by ISO 18562 (all parts).

NOTE 2 This document has been prepared to address the relevant essential principles of safety and performance as indicated in Annex B.

2 Normative reference

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 7396-1:2016, *Medical gas pipeline systems — Part 1: Pipeline systems for compressed medical gases and vacuum*

ISO 14971:2007, *Medical devices — Application of risk management to medical devices*

ISO 18562-1:2017, *Biocompatibility evaluation of breathing gas pathways in healthcare applications — Part 1: Evaluation and testing within a risk management process*

3 Term and definition

For the purposes of this document, the terms and definitions given in ISO 7396-1, ISO 14971, ISO 18562-1 and the following apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

— IEC Electropedia: available at <http://www.electropedia.org/>

— ISO Online browsing platform: available at <http://www.iso.org/obp>

3.1

DIAMETER

AERODYNAMIC DIAMETER

diameter of a sphere of density 1 g/cm³ with the same terminal velocity due to gravitational force in calm air as the particle of interest, regardless of its geometric size, shape and true density, under the prevailing conditions of temperature, pressure and relative humidity

[SOURCE: ISO 7708:1995, 2.2, modified — added “of interest, regardless of its geometric size, shape and true density”]

4 General principles

4.1 TYPE TESTS

The tests described in this document are TYPE TESTS. TYPE TESTS are performed on the final MEDICAL DEVICE, a component of the MEDICAL DEVICE or a representative sample of the MEDICAL DEVICE, part or ACCESSORY being evaluated. If representative samples are used (i.e. manufactured and processed by equivalent methods), consideration should be given to whether or not the differences between the representative sample and the final MEDICAL DEVICE or component could affect the results of the test. Testing of representative samples (manufactured and processed by equivalent methods) instead of the final MEDICAL DEVICE should be supported by a description of any differences between the representative sample and the final MEDICAL DEVICE, and a detailed rationale for why each difference is not expected to impact the BIOCOMPATIBILITY of the final MEDICAL DEVICE.

4.2 General

All GAS PATHWAYS from which the PATIENT inspires gas shall be evaluated using the strategy detailed in ISO 18562-1.

5 * PARTICULATE MATTER emissions

5.1 General

A MEDICAL DEVICE, part or ACCESSORY shall not add to the gas that could be inspired by the PATIENT levels of PARTICULATE MATTER:

- less than or equal to 2,5 µm DIAMETER, in excess of 12 µg/m³;
- less than or equal to 10 µm DIAMETER, in excess of 150 µg/m³.

NOTE 1 The allowable limits are taken from the US EPA 40 § CFR Part 50^[5].

All GAS PATHWAYS from which the PATIENT inspires gas shall be evaluated for PARTICULATE MATTER emissions. The evaluation should use the RISK MANAGEMENT PROCESS to assess if testing is required.

NOTE 2 The evaluation of some components, which are identical in FORMULATION, processing and preparation for use to an existing component of a MEDICAL DEVICE that has been previously tested, might conclude that no further testing is required. Refer to ISO 18562-1:2017, Figure 2.

Evaluation and, if required, testing shall take in to account:

- the EXPECTED SERVICE LIFE;
- the effects of any intended processing or reprocessing;
- the worst-case PATIENT exposure.

The MANUFACTURER shall document this evaluation as well as the criteria for selection of test articles and methodologies, including component parts to be tested, duration of testing in relation to the intended duration of clinical use.

NOTE 3 Some AUTHORITIES HAVING JURISDICTION evaluate these rationales.

If the RISK MANAGEMENT PROCESS determines that testing is required, the testing according to 5.5, 5.6, or 5.7 shall be performed. For testing according to 5.5, use the setup according to either 5.3 or 5.4. The MANUFACTURER may choose the appropriate test method.

Compliance is checked by RISK MANAGEMENT plan and RISK MANAGEMENT FILE.

5.2 Testing methods overview

There is a great variety of components and MEDICAL DEVICES within the scope of this document, and so several different methods are proposed. The MANUFACTURER should select the most appropriate method for their particular application. A simple component such as a connector with minimal area exposed to the PATIENT breathing gas stream is very unlikely to need testing for PARTICULATE MATTER, while a mechanical MEDICAL DEVICE with moving parts such as a ventilator could well require thorough testing.

The simplest method (described in 5.3) is to use a single particle filter to trap everything with a DIAMETER over $0,2\ \mu\text{m}$, and consider the limit to be $12\ \mu\text{g}/\text{m}^3$ for all trapped particles. This is a quick simple test that does not differentiate particle sizes. It may be sufficient for simple MEDICAL DEVICES. It is very difficult to measure very small amounts of PARTICULATE MATTER captured using a barrier filter test method since the mass of the filter is substantially more than that of the PARTICULATE MATTER. The volume of gas used in the test should therefore be large enough to capture a sufficiently large amount of PARTICULATE MATTER to be able to measure it or prove that the total mass of PARTICULATE MATTER is below the allowed amount.

If the MANUFACTURER wishes to test for the different particle sizes, with the different limits as detailed from the US EPA 40 § CFR Part 50^[5], then a full test using inertial particle separators and filters following the general principles described in 40 § CFR Part 50 is required. This is described in more detail in 5.6.

A third alternative is to use a particle counter. The particle count measured by these instruments needs to be converted into an estimate of $\mu\text{g}/\text{m}^3$. A method is suggested in 5.7.


5.3 Single filter test setup

This is a simple method suitable for suspected low levels of PARTICULATE MATTER.

In principle, if sufficiently clean input gas is available, then a single measurement of PARTICULATE MATTER contamination in the output gas stream is sufficient. All of the PARTICULATE MATTER measured is considered to have come from the MEDICAL DEVICE itself as indicated in Figure 1. For a simple, low flow MEDICAL DEVICE, this may be sufficient.

NOTE It is important to ensure that the filter is validated for filtration of particles in airstreams, and that it is suitable for the airflow being used.

The input gas stream may be cleaned by passing all the input air through a $0,2\ \mu\text{m}$ filter before the MEDICAL DEVICE. Then the measuring filter on the output only measures PARTICULATE MATTER that originates from the MEDICAL DEVICE itself.

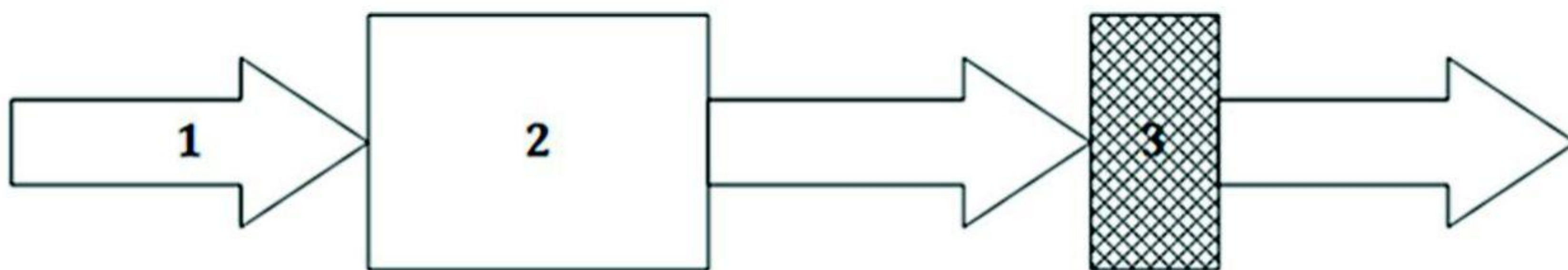

Key


- 1 clean input airstream, filtered if necessary
- 2 one or more MEDICAL DEVICES under test
- 3 $0,2\ \mu\text{m}$ filter

To produce a meaningful result, more than one MEDICAL DEVICE may be required to be placed in series or measured sequentially.

Figure 1 — Example test setup for full flow

If the MEDICAL DEVICE operates with a flowrate, in excess of that which reasonably dimensioned filters can handle, then a different approach may be utilized. For these flowrates, it is not feasible to have the full flow pass through the $0,2\ \mu\text{m}$ filter, so a fractional sampling method is

used as indicated in Figure 2. A subatmospheric pressure (partial vacuum) source may be used to draw the sample volume through the measurement filter

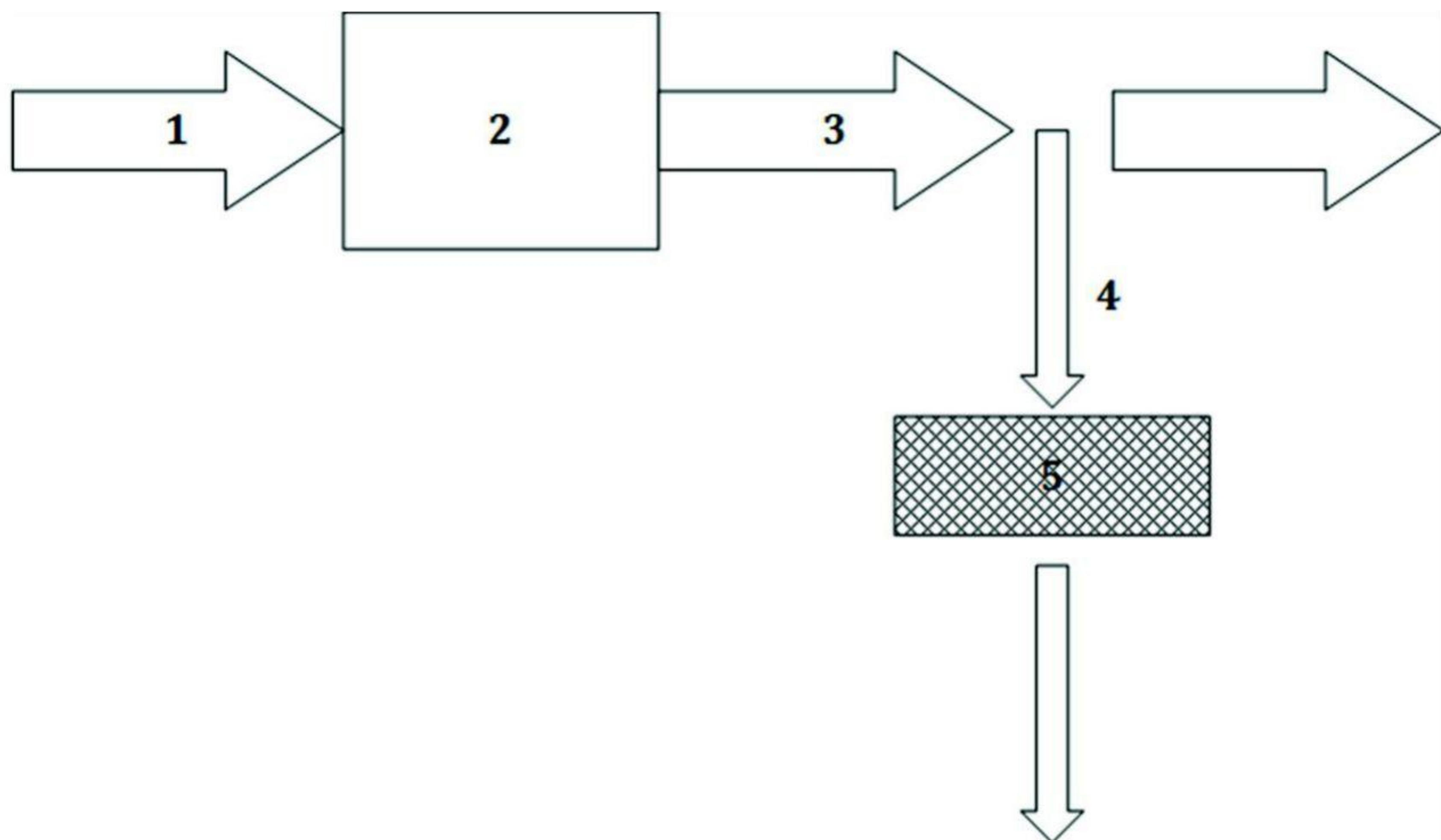

Key


- 1 clean input airstream
- 2 one or more MEDICAL DEVICES under test
- 3 output airstream
- 4 sampled airstream
- 5 0,2 μm filter

To produce a meaningful result, more than one MEDICAL DEVICE may be required to be placed in series or measured sequentially.

Figure 2 — Example single filter test setup for a sampled flow

5.4 Double filter test setup

If sufficiently clean input air is not available for a MEDICAL DEVICE, then a double sampling technique may be used. The principle is to measure the amount of PARTICULATE MATTER in the airstream entering the MEDICAL DEVICE (measurement 1) and simultaneously measure the amount of PARTICULATE MATTER in the airstream leaving the MEDICAL DEVICE (measurement 2), and then subtract the first from the second to get the amount of PARTICULATE MATTER added to the airstream by the MEDICAL DEVICE itself. This method is indicated in Figure 3. A subatmospheric pressure (partial vacuum) source may be used to draw the sample volume through the measurement filters.

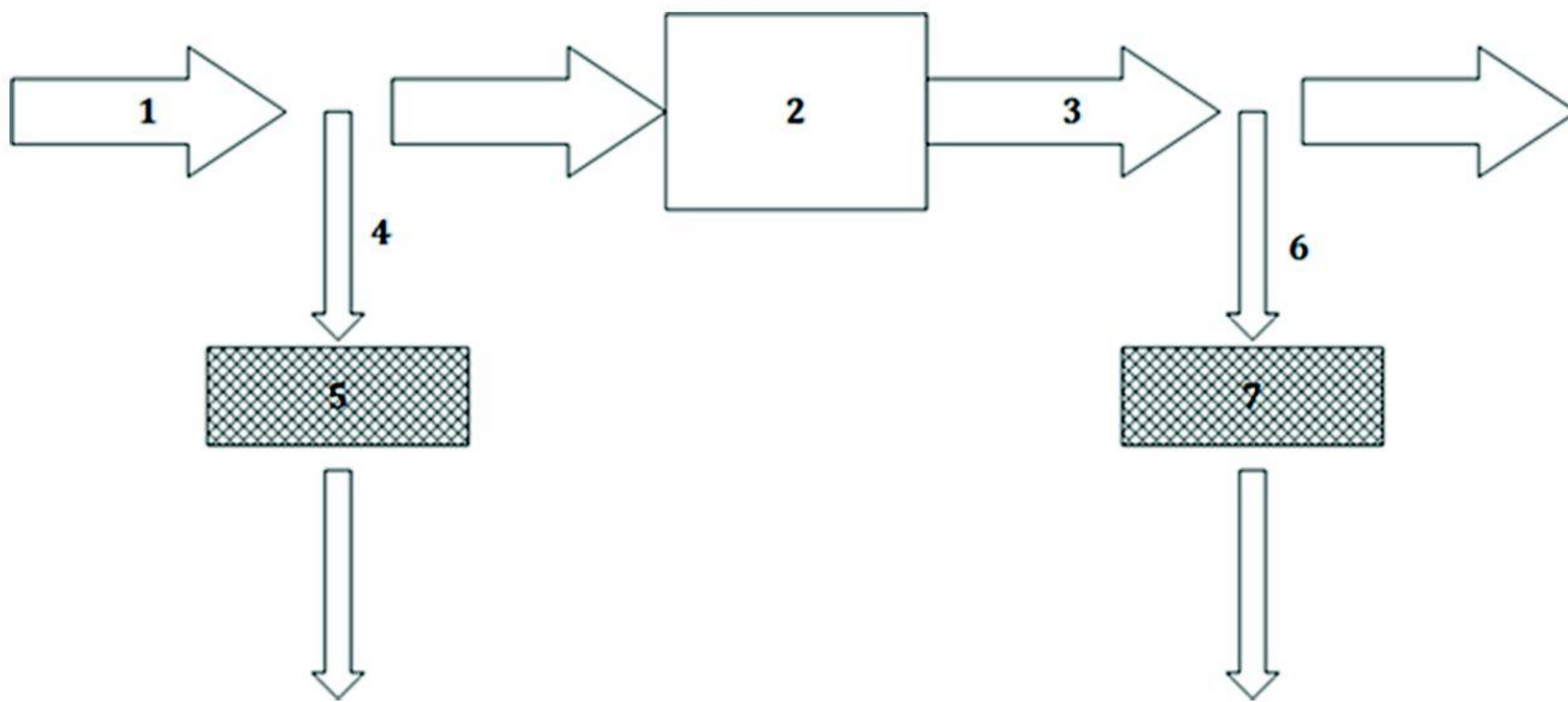

Key

- 1 input airstream
- 2 MEDICAL DEVICE under test
- 3 output airstream
- 4 sampled airstream 1
- 5 0,2 µm filter, measurement 1
- 6 sampled airstream 2
- 7 0,2 µm filter, measurement 2

To produce a meaningful result, more than one MEDICAL DEVICE may be required to be placed in series or measured sequentially.

Figure 3 — Example double filter test setup for a sampled flow

5.5 Test method

Perform filter method PARTICULATE MATTER emission testing as follows.

- a) Choose filters that are suitable for the flowrates passing through them.
- b) Weigh the 0,2 µm rated filter.
- c) Operate the MEDICAL DEVICE in NORMAL USE at the maximum clinically significant flow.
- d) Determine the flow of gas passing through the filter.
- e) For short duration of use MEDICAL DEVICES (e.g. manual emergency resuscitators), the test duration should be the maximum expected duration of use in the clinical application. More than one MEDICAL DEVICE may be needed to create adequate volume to allow reliable measurement. For long duration of use MEDICAL DEVICES, adjust the duration of the test to ensure that the sampling volume through the filter is large enough to permit reliable measurement of 2,5 µg/m³.

NOTE 1 This is approximately 20 % of the 12 µg/m³ limit.

EXAMPLE For a balance that can accurately measure a 20 µg increase in mass of a filter, sample long enough to pass 8 m³ through the filter (8 m³ × 2,5 µg/m³ = 20 µg). This provides an adequate margin to ensure that the result is meaningful.

- f) Re-weigh the filter.
- g) Subtract the mass measured before the test began from the mass after the test is finished to determine the added mass of the particles filtered out of the gas stream by that filter.
- h) Confirm that the level of PARTICULATE MATTER emitted by the MEDICAL DEVICE is less than the limit specified in 5.1.

NOTE 2 Since this method does not differentiate between particle sizes, then the limit of 12 $\mu\text{g}/\text{m}^3$ is used.

NOTE 3 Care needs to be taken to ensure that the moisture content of the filter is the same before and after the test.

NOTE 4 Care needs to be taken that there is no liquid water condensed out on the filter, which might cause loss of effective filtering.

NOTE 5 For the double filter method, care needs to be taken to ensure that the same volume of gas passes through each filter.

5.6 Measuring PARTICULATE MATTER emissions according to particle size

If the MANUFACTURER wishes to test for the different particle sizes, with the different limits as detailed from the US EPA 40 \S CFR Part 50^[5], then a full test using inertial particle separators and filters according to the methods in 40 \S CFR Part 50 is required.

Confirm that the level of PARTICULATE MATTER emitted by the MEDICAL DEVICE is less than the limit specified in 5.1.

5.7 * Measuring PARTICULATE MATTER emissions by particle counter

As an alternative to the filter methods, a calibrated particle counter may be used to measure the PARTICULATE MATTER in the airstream emerging from the MEDICAL DEVICE. The particle counter gives a count of the number of particles of a particular size detected. This number of particles needs to be converted to a mass (in $\mu\text{g}/\text{m}^3$) of air.

NOTE 1 Additional information is found in Annex A.

NOTE 2 Some AUTHORITIES HAVING JURISDICTION require the MANUFACTURER to validate this test method against a filter test method.

Confirm that the level of PARTICULATE MATTER emitted by the MEDICAL DEVICE is less than the limit specified in 5.1.

5.8 MEDICAL DEVICES with time varying emissions

For MEDICAL DEVICES or components whose PARTICULATE MATTER emissions can change over time, it might be necessary to repeat tests for PARTICULATE MATTER after simulating appropriate periods of use, as determined by the RISK MANAGEMENT PROCESS.

Annex A (informative) **Rationale and guidance**

A.1 General guidance

This annex provides rationale for the important requirements of this document and is intended for those who are familiar with the subject of this document but who have not participated in its development. An understanding of the reasons for the main requirements is considered to be essential for its proper application. Furthermore, as clinical practice and technology change, it is believed that rationale for the present requirements will facilitate any revision of this document necessitated by those developments.

The clauses and subclauses in this annex have been so numbered to correspond to the clauses and subclauses in this document to which they refer. The numbering is, therefore, not consecutive.

A.2 Rationale for particular clauses and subclauses

Clause 5 — PARTICULATE MATTER emissions

Penetration and deposition of PARTICULATE MATTER in the human respiratory tract are complicated subjects.

The following is adapted from an extract from Reference [6].

For better understanding, a schematic representation of the respiratory system is presented in Figure A.1. It shows the different regions of the human respiratory system, namely, nasopharyngeal (or extrathoracic), tracheobronchial and alveolar.


PARTICULATE MATTER small enough to stay airborne can be inhaled through the nose (nasal route) or the mouth (oral route). The probability of inhalation depends on the PARTICULATE MATTER DIAMETER (particle AERODYNAMIC DIAMETER), air movement around the body, and breathing rate. The inhaled PARTICULATE MATTER can then either be deposited or exhaled again, depending on a whole range of physiological and PARTICULATE MATTER-related factors. The five deposition mechanisms are sedimentation, inertial impaction, diffusion (significant only for very small PARTICULATE MATTER $<0,5\ \mu\text{m}$), interception and electrostatic deposition. Sedimentation and impaction are the most important mechanisms in relation to inhaled airborne dust, and these PROCESSES are governed by PARTICULATE MATTER DIAMETER. There are significant differences between individuals in the amount deposited in different regions.

The US EPA limits for PARTICULATE MATTER are specified as overall mass of particulates (of specified particle size range) per volume of air, irrespective of the material composition. At present, the scientific literature on HAZARDS of PARTICULATE MATTER supports this approach. Studies that have investigated the properties of inhaled particulates and their relation to health outcomes have identified particle size as strongly correlated with HAZARD to health^{[7][8]}.

The largest inhaled PARTICULATE MATTER, with DIAMETER greater than about $30\ \mu\text{m}$, is deposited in the air passages between the point of entry at the lips or nares and the larynx. During nasal breathing, PARTICULATE MATTER is deposited in the nose by filtration by the nasal hairs and impaction where the airflow changes direction. Retention after deposition is helped

by mucus, which lines the nose. In most cases, the nasal route is a more efficient PARTICULATE MATTER filter than the oral route, especially at low and moderate flowrates. Thus, people who normally breathe part or all of the time through the mouth can be expected to have more PARTICULATE MATTER reach the lung and deposit there than those who breathe entirely through the nose. During exertion, the flow resistance of the nasal passages causes a shift to mouth breathing in almost all people. Other factors influencing the deposition and retention of PARTICULATE MATTER include cigarette smoking and lung disease.

Of the PARTICULATE MATTER, which fails to deposit in the nasopharyngeal region, the larger sizes deposit in the tracheobronchial airway region and can later be eliminated by mucociliary clearance or, if soluble, can enter the body by dissolution. The smaller sizes can penetrate to the alveolar region (see Figure A.1), the region where inhaled gases can be absorbed by the blood. In DIAMETER terms, only about 1 % of 10 μm particles get as far as the alveolar region, so 10 μm usually is considered the practical upper size limit for penetration to this region. Maximum deposition in the alveolar region occurs for particles of approximately 2 μm DIAMETER. Most PARTICULATE MATTER larger than this is deposited further up in the lung. For smaller PARTICULATE MATTER, most deposition mechanisms become less efficient, so deposition is less for PARTICULATE MATTER smaller than 2 μm until it is only about 10 % to 15 % at about 0,5 μm . Most of this PARTICULATE MATTER is exhaled again without being deposited. For still smaller PARTICULATE MATTER, diffusion is an effective mechanism and deposition probability is higher. Deposition is therefore at a minimum at about a DIAMETER of 0,5 μm

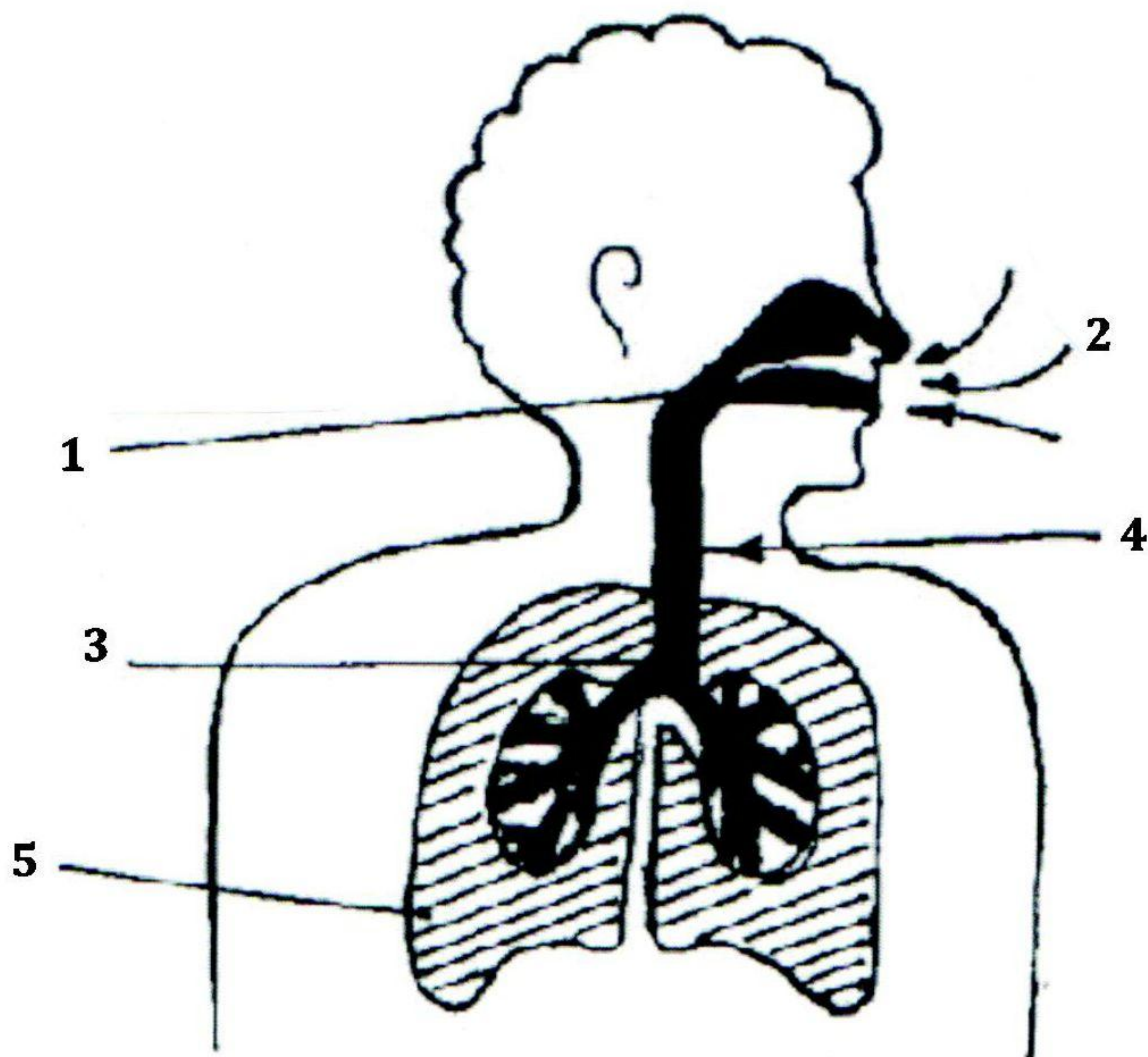

Key

- 1 nasopharyngeal (extrathoracic) region
- 2 inhaled air
- 3 tracheobronchial region
- 4 larynx
- 5 alveolar region

Figure A.1 — Schematic representation of the human respiratory tract

The smaller the PARTICULATE MATTER, the deeper into the lungs it can penetrate and the longer it takes for the body to eliminate it. Originally, the main health concerns with regard to PARTICULATE MATTER were focussed on respiratory health, but now there is emerging evidence of effects on the cardiovascular system as well.

Refer to 40 § CFR Part 50 from the Environmental Protection Agency of the USA for a comprehensive rationale about the health effects of PARTICULATE MATTER and acceptable levels^[5].

Elevated levels of PARTICULATE MATTER exposure have been associated with the declines in lung function and with increases in respiratory system distress such as cough, shortness of breath and asthma attack^[9]. Smaller PARTICULATE MATTER poses greater RISK to health than larger PARTICULATE MATTER because smaller PARTICULATE MATTER is more toxic and is

breathed more deeply into the lungs. Smaller sized PARTICULATE MATTER is retained in the alveolar region and can penetrate even deeper into interstitial sites^{[10][11]}.

The following is provided for information only and to put the particle concentrations mentioned in this document into perspective. Table A.1 is taken from ISO 16000-37:—, Annex B. This table is an informative listing of empirical values obtained for concentration ranges of the fractions PM₁₀, PM_{2,5} and ultrafine particles through indoor air measurements of residential premises in Germany.

Table A.1 — Empirical values for particle concentration ranges of the fractions PM₁₀ and PM_{2,5}

Indoor situation	Measured particle/ fraction	Empirical values of typical concentration ranges µg/m³	Concentration depends in particular on
I Presence and general activities of persons			
Dwellings	PM ₁₀	10 to 80	Number and activity
	PM _{2,5}	10 to 40	
Schools, day nurseries	PM ₁₀	40 to 150	
	PM _{2,5}	10 to 40	
Offices	PM ₁₀	20 to 60	
	PM _{2,5}	10 to 40	
II Specific user activities			
Smoking	PM ₁₀	50 to 500	Number/quantity
	PM _{2,5}	20 to 100	
Using a vacuum cleaner	PM ₁₀	30 to 150	Degree of pollution, filtration performance
	PM _{2,5}	10 to 40	
Cooking/preparing hot water	PM ₁₀	40 to 100	Duration and intensity
Stove/fireplace	PM ₁₀	40 to 200	Fireplace/stove construction, heating material, chimney

Subclause 5.7 — Measuring PARTICULATE MATTER emissions by particle counter

Refer to the instructions for use of the particle counter to clarify exactly what the output reading of the particle counter means and in what units it is presented.

To convert from “number of particles per cubic metre” to a “mass per cubic metre”, the average density of the particles needs to be determined. For most MEDICAL DEVICES, this can be estimated with sufficient accuracy for the purposes of this document by considering the materials of manufacture of the MEDICAL DEVICE. For example, if the MEDICAL DEVICE is of mostly plastic construction, then an average density of the type of plastics used in the construction can be used, since emitted particles will be most likely to arise from plastic base materials. In these calculations, the shape of the particle is assumed to be spherical, in keeping with the definition of DIAMETER.

EXAMPLE A typical density of a polymer used in breathing circuits is 0,9 g/cm³. So a particle with a DIAMETER of 1 µm has a volume of $[(4/3) \cdot \pi \cdot r^3] = 0,52 \mu\text{m}^3, = 0,52 \times 10^{-12} \text{ cm}^3$ and would therefore weigh $0,52 \times 0,9 \times 10^{-12} \text{ g} = 0,47 \times 10^{-12} \text{ g} = 0,47 \times 10^{-6} \mu\text{g}$. So each particle weighs this much. If the particle counter gives a reading of 10⁶ particles/m³ (DIAMETER of 1 µm), then using the mass derived above, the mass of 1 µm particles/m³ is $0,47 \times 10^{-6} \times 10^6 \mu\text{g} = 0,47 \mu\text{g}$.

Other MEDICAL DEVICES are made of other materials, for example, aluminium for a turbine motor housing. The particles likely to arise from such a MEDICAL DEVICE would then have the density of aluminium, and this density (approximately 2,7 g/cm³) is used in the calculation.

If the source or composition of the particles emitted is not known, then this density to mass conversion is not easily accomplished. In this case, the worst-case density of materials in the MEDICAL DEVICE is used. To simplify the calculations, it is possible to assume that the particles emitted are all of the densest material likely to arise from the MEDICAL DEVICE.

Annex B
(informative)
Reference to the essential principles

This document has been prepared to support the essential principles of safety and performance of GAS PATHWAYS as components of MEDICAL DEVICES according to ISO 16142-1^[3]. This document is intended to be acceptable for conformity assessment purposes.

Compliance with this document provides one means of demonstrating conformance with the specific essential principles of ISO 16142-1. Other means are possible. Table B.1 maps the clauses and subclauses of this document with the essential principles of ISO 16142-1.

Table B.1 — Correspondence between this document and the essential principles

Essential principle of ISO 16142-1	Corresponding clause(s)/ subclause(s) of this document	Qualifying remarks/notes
8.1 a)	Clause 4, Clause 5	Only the part relating to toxicity is addressed.
8.1 b)	Clause 4, Clause 5	
8.2	Clause 4, Clause 5	
8.4	Clause 4, Clause 5	
8.5	Clause 4, Clause 5	Only the part relating to egress of substances from the MEDICAL DEVICE is addressed.

Annex C
(informative)
Terminology — Alphabetized index of defined terms

NOTE The ISO Online Browsing Platform (OBP) and the IEC Electropedia provide access to many of these terms and definitions.

Term	Source
AERODYNAMIC DIAMETER	3.1
ACCESSORY	ISO 18562-1:2017, 3.1
AUTHORITY HAVING JURISDICTION	ISO 16142-1:2016, 3.1
BIOCOMPATIBILITY	ISO 18562-1:2017, 3.2
DIAMETER	3.1
EXPECTED SERVICE LIFE	ISO 18562-1:2017, 3.3
FORMULATION	ISO 18562-1:2017, 3.4
GAS PATHWAY	ISO 18562-1:2017, 3.5
HAZARD	ISO 14971:2007, 2.2
MANUFACTURER	ISO 14971:2007, 2.8
MEDICAL DEVICE	ISO 18562-1:2017, 3.7
MEDICAL GAS PIPELINE SYSTEM	ISO 7396-1:2016, 3.36
NORMAL USE	ISO 18562-1:2017, 3.9
PARTICULATE MATTER	ISO 18562-1:2017, 3.10
PATIENT	ISO 18562-1:2017, 3.11
PROCESS	ISO 14971:2007, 2.13
RISK	ISO 14971:2007, 2.16
RISK MANAGEMENT	ISO 14971:2007, 2.22
RISK MANAGEMENT FILE	ISO 14971:2007, 2.23
TYPE TEST	ISO 18562-1:2017, 3.15

Bibliography

- [1] ISO 7708:1995, *Air quality — Particle size fraction definitions for health-related sampling*
- [2] ISO 10993 (all parts), *Biological evaluation of medical devices*
- [3] ISO 16142-1:2016, *Medical devices — Recognized essential principles of safety and performance of medical devices — Part 1: General essential principles and additional specific essential principles for all non-IVD medical devices and guidance on the selection of standards*
- [4] ISO 16000-37:—¹, *Indoor air — Part 37: Strategies for the measurement of PM 2.5*
- [5] US Code of Federal Regulations, 40 § CFR Part 50, National primary and secondary ambient air quality standards. Available from: http://www.ecfr.gov/cgi-bin/text-idx?tpl=/ecfrbrowse/Title40/40cfr50_main_02.tpl
- [6] World Health Organization (WHO), *Indoor air quality: Organic pollutants. EURO Reports and Studies*, No. 111. Copenhagen
- [7] HARRISON R.M. Particulate matter in the atmosphere: which particle properties are important for its effects on health. *Sci. Total Environ.* Apr 2000, **249**(1) pp. 85–101
- [8] SCHWARZE P.E. et al. Particulate matter properties and health effects: consistency of epidemiological and toxicological studies. *Hum. Exp. Toxicol.* Oct 2006, **25**(10) pp. 555–579
- [9] WHO/SDE/OEH, *Hazard Prevention and Control in the Work Environment: Airborne Dust*, 99.14. Available from: http://www.who.int/occupational_health/publications/en/oehairbornedust3.pdf
- [10] DOCKERY D.W., POPE C.A., XU X.P., SPENGLER J.D., WARE J.H., FAY M.E. et al. An association between air-pollution and mortality in 6 United-States cities. *N. Engl. J. Med.* 1993, pp. 329
- [11] OBERDORSTER G., OBERDORSTER E., OBERDORSTER J. Nanotoxicology: An emerging discipline evolving from studies of ultrafine particles. *Environ. Health Perspect.* 2005, pp. 113

¹ Under preparation. Stage at the time of publication: ISO/AWI 16000-37:2017.

Informasi pendukung terkait perumus standar

[1] Komite Teknis

Komite Teknis 11-03 Alat Kesehatan Elektromedik

[2] Susunan Keanggotaan Komite Teknis

Ketua	: Marlina Harahap	— Balai Pengamanan Fasilitas Kesehatan Jakarta
Wakil Ketua	: Hendrana Tjahjadi	— Asosiasi Perusahaan Laboratorium Pengujian dan Kalibrasi Fasilitas Kesehatan
Sekretaris	: Amjad Tri Puspitasari	— Badan Standardisasi Nasional
Anggota	: 1. Jojo Simanjuntak	— Kementerian Kesehatan
	2. Rakhmat Sauma	— Gakeslab
	3. Chasri Idham	— Asosiasi Produsen Alat Kesehatan Indonesia
	4. Agus Komarudin	— Ikatan Elektromedis Indonesia
	5. Arif Jauhari	— Politeknik Kesehatan Kementerian Kesehatan Jakarta II
	6. Ahmad Bilal	— Perhimpunan Rumah Sakit Seluruh Indonesia

[3] Konseptor

Gugus Kerja Komtek 11-03

[4] Sekretariat Pengelola Komite Teknis

Direktorat Pengembangan Standar Agro, Kimia, Kesehatan dan Halal
Badan Standardisasi Nasional